



Fidelis

Secondary Standard Dose Calibrator
Software Installation Notes
Issue 1

Contents

Revision History	2
Installation Notes	3
Microsoft .Net Framework.....	3
Fidelis Installation.....	4

Revision History

Issue	Date	Changed By
1	April 2007	First Edition

Installation Notes

The Installation CD can be used to install the Fidelis application software.

It can also be used to reinstall the calibration factor database.

After installation, the following folders and files will be installed to your C:\ drive:


Folder	File	
C:\fidelis	Fidelis.exe	Fidelis executable
	Release.txt	Release notes
C:\fidelis\database	fidelis_database_main.csv	main NPL verified calibration factor database
	fidelis_database_tests.csv	daily test results and pass conditions database
C:\fidelis\manual	Fidelis User Manual	Fidelis user manual in PDF format
	Issue3.pdf	
C:\fidelis\drivers\usb	CDM_Setup.exe	USB driver files required for Win2000 & XP

Microsoft .Net Framework


The Fidelis Application requires Microsoft .NET Framework Version 2.0 (x86) to be installed. This is already installed when the system is delivered, but if you need to re-install the system software for any reason, you can find this package in the drivers\DotNet2 folder on the Installation CD.

Fidelis Installation


Click Setup.exe file found in the root folder of the Fidelis Installation CD.

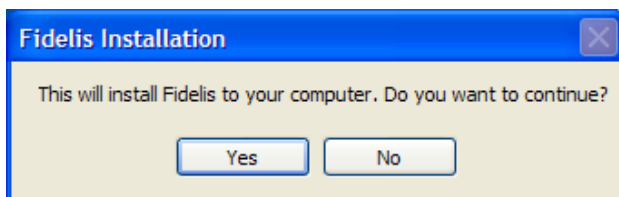


(1) Click Yes to continue.

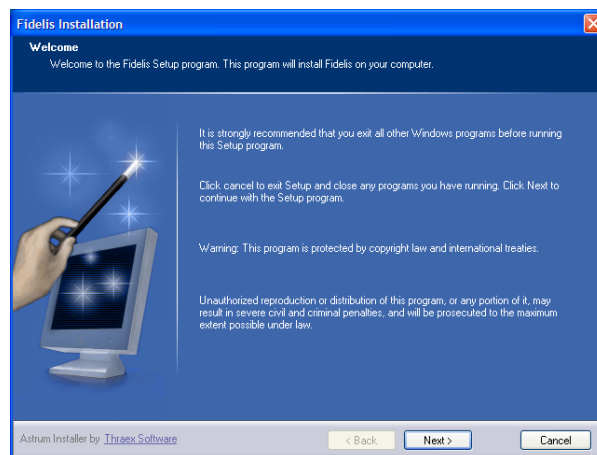


(2) Click Next

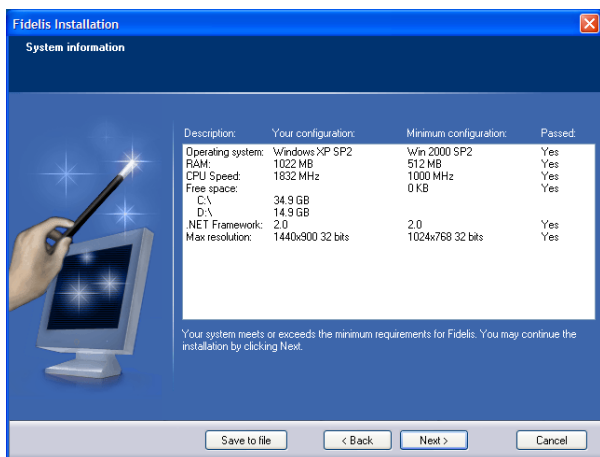

(3) Click Next to confirm that your computer meets the minimum recommended specification to run the Fidelis software.

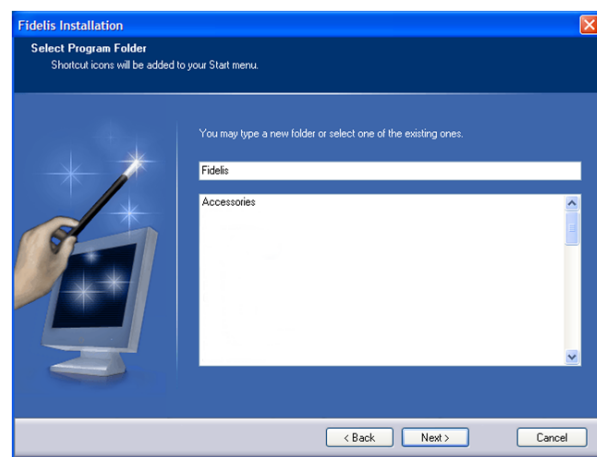


(4) Click next to confirm the Start Menu Program Folder name.

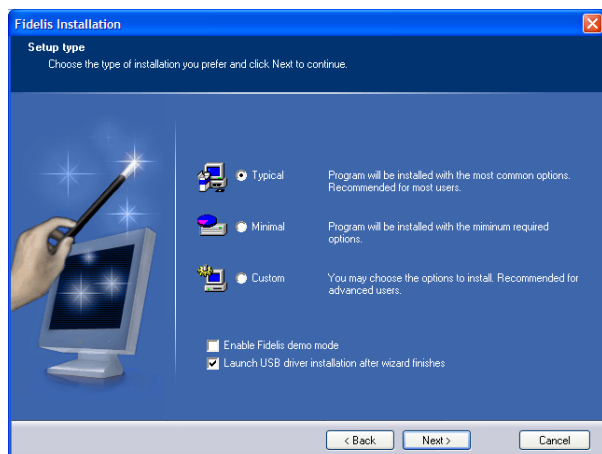


- (5) Select Typical to install both the application Fidelis software and the main calibration factor database. Select Custom if you want to install just the database.

Selecting the “Enable Fidelis demo mode” option will instruct the installer to place the Fidelis software in Demonstration Mode. This allows the software to be used with having a chamber connected.

Uncheck the Launch USB driver install to prevent this running after you finish the wizard.

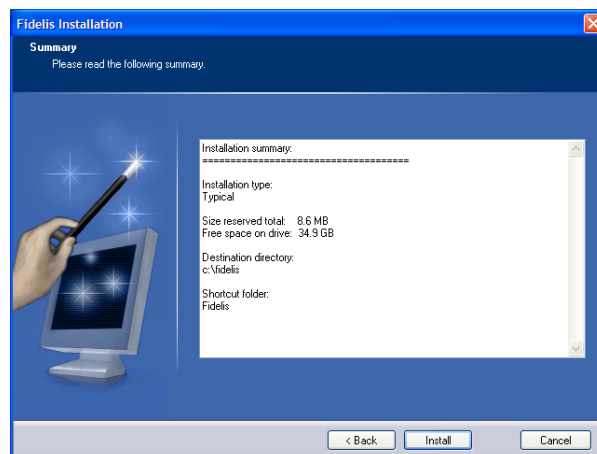


- (6) Click install to start the installation

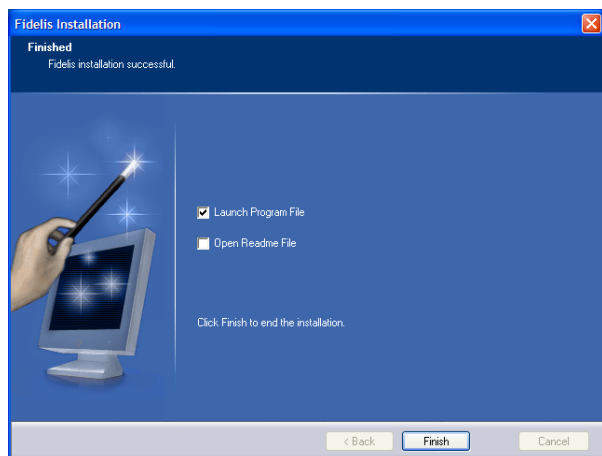

- (7) Click Finish to close the installation wizard and (optionally) start installation of the USB driver files.

Wait until the next dialog appears indicating that the FTDI CDM USB drivers have been installed.

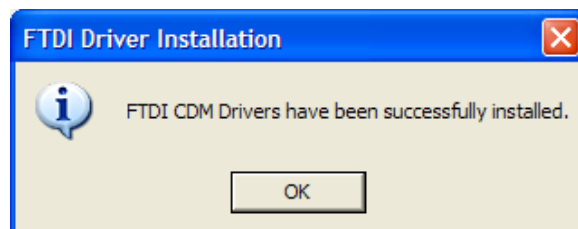

- (8) Click OK to confirm the USB driver installation task.